


MMSI



McARTHUR
MEDICAL SALES INC.

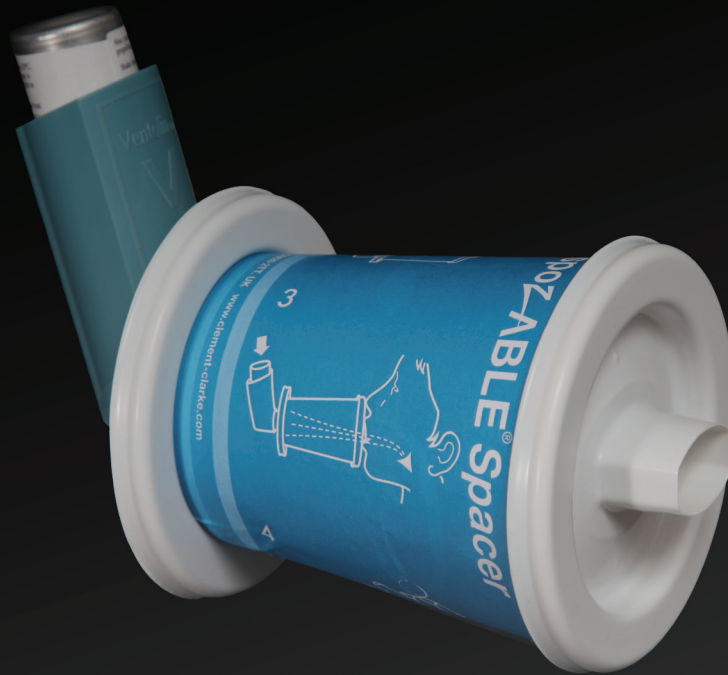

DispozABLE Spacer

Technical specifications

Tradition and Innovation - Since 1858, visionary thinking and a fascination with technology have guided us to develop innovative products of outstanding reliability. Anticipating trends to improve the quality of life.

 **CLEMENT CLARKE**
INTERNATIONAL

DispozABLE Spacer

Technical specifications

DispozABLE™ Spacer

DispozABLE Spacer

| | |
|--------------------|---------------|
| Dimensions | ø84mm x 107mm |
| Weight | 15g |
| Volume | 227ml |
| Resistance to flow | Free flow |

Performance Data

Fine particle dose (salbutamol):

Delay 0 secs = 54.5mcg (SD 4.8)

Delay 1 secs = 37.7mcg (SD 7.6)

| Product | Part Code |
|-------------------------------|-----------|
| DispozABLE Spacer (box of 10) | 3609550 |

Product information

DispozABLE Spacer has been designed to meet a specific market need. It is a low cost, fast response spacer for delivering medication promptly. It is designed to be used straight from the box without pre-washing in an emergency situation. For example in an ambulance or school. DispozABLE Spacer is to be used as a single use spacer. This simple intuitive design allows non-professionals to assemble it quickly. Clear instructions are included on the spacer, adding to the element of ease.

- Simple to use
- Effective
- For use in an emergency situation: ambulance or school
- Easy to assemble
- No mask required
- Can be used with tidal breathing technique, assisted by Flo-Tone.

