

Vitalograph Respiratory Test Disposables



the difference is clear








Respiratory Test Disposables from Vitalograph

Vitalograph[®]
Your respiratory partner

BVF™ Bacterial/Viral Filters


- Single use filters for effective removal of bacteria, viruses and expectorated matter
- Suitable for inspiratory and expiratory spirometry testing
- Protects the subject against cross-infection
- Protects the spirometer against cross-contamination
- Validated 99.99% cross-contamination efficiency for both bacteria and viruses
- Very low expiratory flow impedance
- Low effective dead space of 65ml as required for lung volume and DLCO testing
- Efficient at high flow rates
- Attractively priced



Code:	Colour:	To fit spirometers from:
28350	Blue  (Small Mouthpiece)	Vitalograph, MIR, Micro Medical and most office spirometers
28362	White 	Most Lab PF Equipment, including Sensormedics, PK Morgan, Ganshorn (new)
28363	Clear 	Jaeger/VIASYS/Carefusion PF Lab Equipment
28365	Teal 	MedGraphics PF Lab Equipment
28375	Grey 	Chest, Fukuda, Minato
28368	Green 	Ferraris/nSpire/PDS/KoKo
28401	Translucent  (Fits Syringe)	Most office spirometers, including Vitalograph, Medisoft, Cosmed, MIR, Micro Medical, ZAN, H.R and other applications

Unlike barrier filters which trap expectorated matter whilst allowing viruses and bacteria to pass through, Vitalograph filters use an electrostatically charged material that traps expectorated matter plus bacteria and viruses for highly effective protection against cross-contamination.

In addition to the hygiene benefits of using bacterial viral filters they also support productivity in pulmonary function testing by reducing the time spent cleaning and decontaminating test equipment.

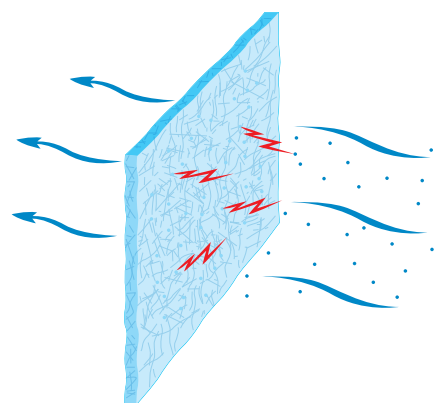

Pulmonary Function Test Kits

- Complete pulmonary function test kits for use with laboratory lung function analysers and office spirometers
- Designed to provide high performing protection against cross-infection
- Ergonomically designed for comfort and ease of use
- Kits can include BVF, optional mouthpiece, optional elbow and noseclip
- Attractively priced
- Completely disposable

All PFT Kits can be supplied with any Vitalograph BVF or mouthpiece, even if not listed overleaf.



Silicone bite-on mouthpiece may be single use or cleaned and disinfected (user's choice)




Synthetic bite-on and plastic mouthpiece, single-use, disposable

Your choice of optional bite-on and plastic mouthpieces

- High quality mouthpieces designed for use with any size filter
- Ergonomically designed for a comfortable fit
- Attractively priced
- Completely disposable or washable

SafeTway® Mouthpieces

- Protects the subject from cross-infection
- Suitable for expiratory manoeuvres with peak flow meters and spirometers
- One-way valve prevents inhalation from spirometer

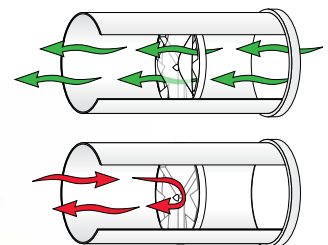

Technical Specification: Ordering Info:

Product Name: B/V Filter

Model No: 2820

Efficiency: 99.99% @ 55 and 750 L/min

Validation using Bacillus Subtillis var niger mean dia.: 0.7 to 1 µm

Resistance to Airflow @ 750L/min (BVF alone): 0.08kPa/L/s (0.8 cmH₂O)

Resistance to Airflow @ 750L/min (BVF + Fleisch): 0.12kPa/L/s (1.2 cmH₂O)

Effective Deadspace (excluding Tapers): 65mL

Material: Electrostatic nanofibre technology filter medium

All Respiratory Test Disposables

Designed & Manufactured Under:
ISO 13485:2003, FDA 21CFR820, CMDR

Medical Safety Standard:
Medical Devices Directive 93/42/EEC
(as amended)

BVF Bacterial/Viral Filters:

28350	2820 Blue BVF fits Office Spirometers (50)
28362	2820 White BVF fits PFT Lab Systems (50)
28363	2820 Clear BVF fits Jaeger (50)
28365	2820 Teal BVF fits Med Graphics (50)
28368	2820 Green BVF fits PDS/Koko (50)
28375	2820 Grey BVF fits Fukuda/Chest/Minato (50)
28401	2820 Translucent BVF fits Cosmed/Medisoft (50)

SafeTway® Mouthpieces:

20242	2024 SafeTway mouthpieces (200)
20980	2024 SafeTway mini mouthpieces (50)

Bite-On Mouthpieces:

20310	2020 Silicone Bite-On Mouthpieces (50)
20989	2020 Disposable Oval Mouthpieces (50)
20305	2020 Disposable Elbow (50)

Noseclips:

20303	2030 Disposable Noseclips (200)
28370	Blue BVF + Noseclip (50)
28371	Clear BVF + Noseclip (50)
28372	White BVF + Noseclip (50)

Some of the available Pulmonary Function Test Kits:

28381	PFT Kit White BVF Silicone (50)
28382	PFT Kit Clear BVF Silicone (50)
28385	PFT Kit Clear BVF Oval Elbow (50)
28400	PFT Kit Translucent BVF Oval (50)
28421	PFT Kit Translucent BVF Silicone (50)
28425	PFT Kit Translucent BVF Disposable (50)

Vitalograph is a world leading provider of outstanding quality cardio- respiratory diagnostic devices, clinical trial services and medical equipment servicing. With a pioneering heritage of excellence spanning half a century Vitalograph continues to make valuable contributions to effective medical care and enhanced quality of life.



1846 5th Concession W., P.O. Box 7, Rockton, Ontario LOR 1X0
Tel.: (519) 622-4030 • Toll Free: 1-800-996-6674 • Fax: (519) 622 - 1142
Email: mmsi@mcarthurmedical.com • Web: mcarthurmedical.com



UK Sales

Vitalograph Ltd.
Maids Moreton Buckingham
MK18 1SW England
Tel: 01280 827110
Fax: 01280 823302
e-mail: sales@vitalograph.co.uk
www.vitalograph.co.uk

Export Sales

Vitalograph Ltd.
Maids Moreton Buckingham
MK18 1SW England
Tel: +44 1280 827120
Fax: +44 1280 823302
e-mail: sales@vitalograph.co.uk
www.vitalograph.eu

Vitalograph GmbH

Rellinger Strasse 64a
20257 Hamburg
Germany
Tel: (040) 547391-0
Fax: (040) 547391-40
e-mail: info@vitalograph.de
www.vitalograph.de

Vitalograph Inc.

13310 West 99th Street
Lenexa
Kansas 66215 USA
Tel: (913) 888 4221
Fax: (913) 888 4259
e-mail: vitcs@vitalograph.com
www.vitalograph.com

Vitalograph (Ireland) Ltd.

Gort Road Business Park
Ennis
Co Clare Ireland
Tel: (065) 6864100
Fax: (065) 6829289
e-mail: sales@vitalograph.ie
www.vitalograph.ie

