


# BlenderBuddy 2

SENSING

ANALYSIS

DELIVERY

*"A compact and clean addition  
to your blender"*

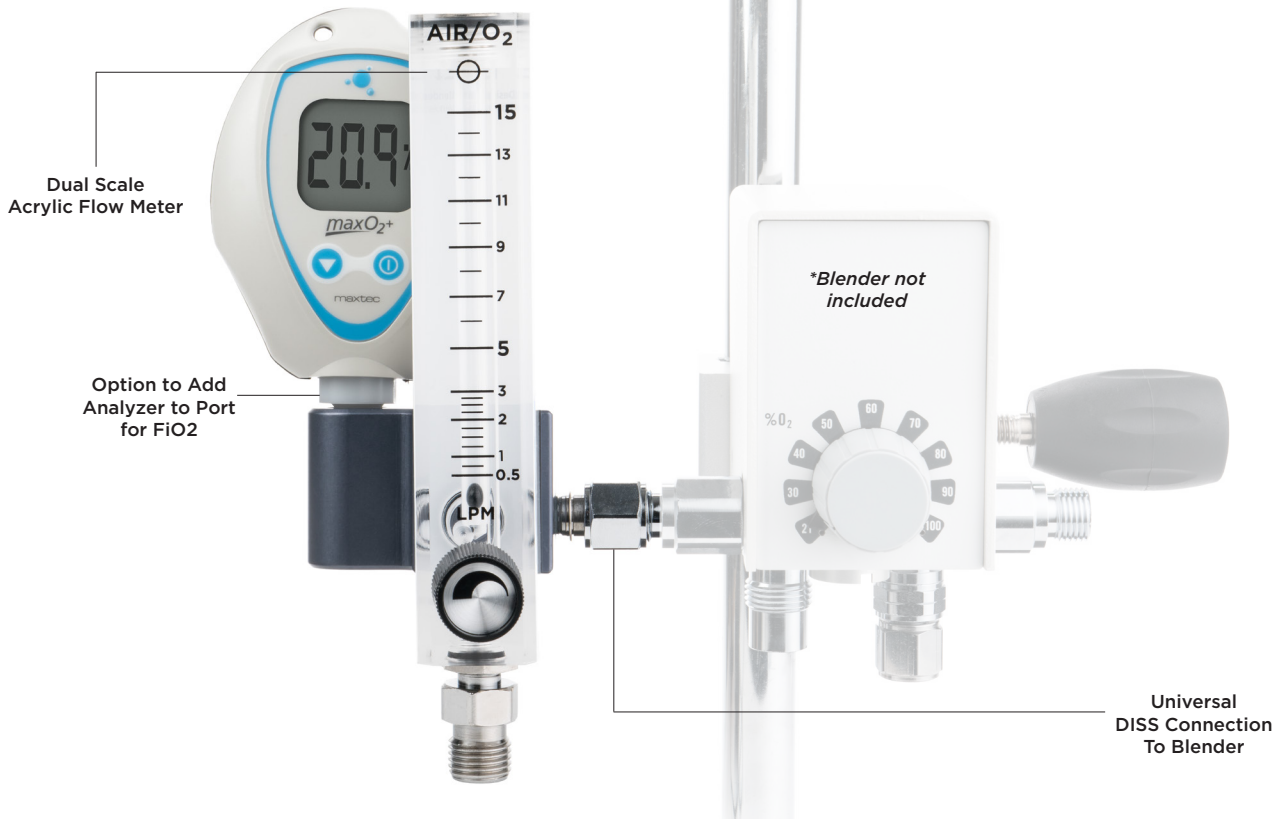

## BLENDERBUDDY 2 OVERVIEW

The **BlenderBuddy 2** equips your blender with the Maxtec **DFB** (*designed for blenders*) dual-scale flow meter and also has a built-in oxygen port, which allows gas to be analyzed safely and away from the patient to avoid cross-contamination. It securely attaches to your blender using the DISS outlet connection. The BlenderBuddy 2 also has a micro bleed, which reduces the amount of wasted gas by more than 15 times previous versions. This is all backed by a **Maxtec three-year warranty**.

## PART NUMBERS

### BlenderBuddy 2, DFB:

0-3 LPM (*Low Flow*) . . . . . R233P01-001  
0-15 LPM (*Low Flow*) . . . . . R233P01-002  
0-30 LPM (*Low Flow*) . . . . . R233P01-003  
0-70 LPM (*High Flow*) . . . . . R233P01-004

### BlenderBuddy 2, DFB, With MaxO2+A:

0-3 LPM (*Low Flow*) . . . . . R233P04-001  
0-15 LPM (*Low Flow*) . . . . . R233P04-002  
0-30 LPM (*Low Flow*) . . . . . R233P04-003  
0-70 LPM (*High Flow*) . . . . . R233P04-004

---

## Specifications

Flow Meter Accuracy	+/- 10% of indicated value or 0.5 lpm, whichever is greater with 50 psig inlet pressure
Device Dimensions (L x W x H)	3.5" x 3.5" x 7"
Device Weight	0.88 lbs (0.4 kg)
Package Dimensions (L x W x H)	8" x 8" x 6"
Package Weight	1.8 lbs (0.8 kg)
Total System Response Time	≤ 30 seconds
Sensor Bleed Flow	≤ 0.2 LPM
Latex Free	Yes
PVC Free	Yes
RoHS Compliant	Yes
Warranty	36 months ( <i>excludes analyzer</i> )



**1846 5th Concession W., P.O. Box 7, Rockton, Ontario LOR 1X0**  
**Tel.: (519) 622-4030 • Toll Free: 1-800-996-6674 • Fax: (519) 622 - 1142**  
**Email: [mmsi@mcarthurmedical.com](mailto:mmsi@mcarthurmedical.com) • Web: [mcarthurmedical.com](http://mcarthurmedical.com)**